



EASYCAP

EEG Recording Caps and Related Products

EASYCAP GmbH

Am Anger 5

82237 Würthsee

GERMANY


Tel +49 8153 88702-00

www.easycap.de

info@easycap.de

32Ch Standard BrainCap with Multitrodes

Electrode Layout and Channel Assignment

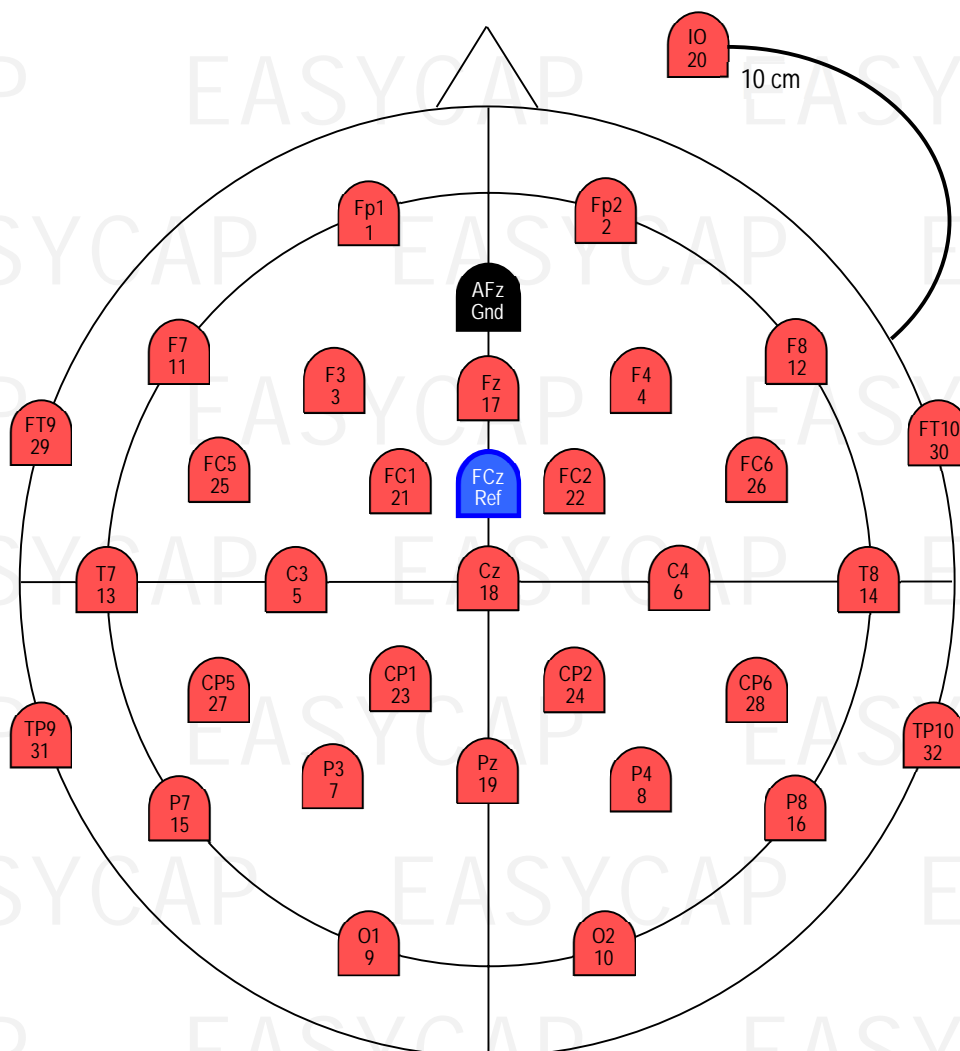

Details

Ordering Information

For ordering please give Article Number, Cap Cut, and Size (e.g. *BC-32, C-Cut, 56*):

- Article Number: *BC-32*
- Cap Cut: *C-Cut or A-Cut*
- Size (given in cm head circumference):
 - Adult caps: *54, 56, 58, 60, 62, 64* (average male: 58, average female: 56)
 - Children caps: *50* (3-4 years), *52* (5-10 years), *54* (11-14 years)
 - Infant caps: *34, 36* (newborn), *38, 40* (3 months), *42, 44* (7 month), *46, 48* (2 years)

The catalogue-number comprises the cap as described. For further information about accessories or consumables, please visit our website or contact our local distributor.

Cap

Standard: Subtemporal Cap with integrated chin belt, white.

Sizes 52 – 60 made from High Precision Fabric, Sizes 50 and smaller made from High Comfort Fabric

Options: C-Cut or A-Cut, Size. For further variations, contact us.

Electrodes

All electrodes are Multitrodes with sintered Ag/AgCl sensors. They are buttoned directly into the cap (total height 3,5 mm) or can be attached to the skin with washers (= double-sided adhesive rings).

All electrodes are number-labelled at the sensor end.

The cable colours correspond to the above figure.

The cables are attached to the cap with nylon threads. Two cable trees leave the cap plait-like behind the ears, pointing downwards, and meet at the connector-box. The length of cable trees is approx. 120 cm.

Termination

The cable trees meet at a Connector box. From here the caps are connected to BrainAmp with 30 cm-flat-ribbon-cables. These flat ribbon cables come with the BrainAmps. They can be re-ordered from BrainProducts (Cat-No. BP-02400-NN) or from EasyCap (Cat.-No. E80).

Theta/Phi-Coordinates

Please find a table with Theta/Phi-Coordinates of all electrode sites at the end of this file.

Table of Coordinates for BC-32

Channel-number	Name	Theta	Phi
1	Fp1	-90	-72
2	Fp2	90	72
3	F3	-60	-51
4	F4	60	51
5	C3	-45	0
6	C4	45	0
7	P3	-60	51
8	P4	60	-51
9	O1	-90	72
10	O2	90	-72
11	F7	-90	-36
12	F8	90	36
13	T7	-90	0
14	T8	90	0
15	P7	-90	36
16	P8	90	-36
17	Fz	45	90
18	Cz	0	0
19	Pz	45	-90
20	IO	---	---
21	FC1	-31	-46
22	FC2	31	46
23	CP1	-31	46
24	CP2	31	-46
25	FC5	-69	-21
26	FC6	69	21
27	CP5	-69	21
28	CP6	69	-21
29	FT9	-113	-18
30	FT10	113	18
31	TP9	-113	18
32	TP10	113	-18
Gnd	Afz	67	90
Ref	FCz	23	90

These values are standardized to a Theta of 90° for the plane through Fpz, T7, T8, Oz.

The signs follow this convention:

