


## 32Ch Standard BrainCap for actiCHamp Plus with Multitrodes

### Electrode Layout and Channel Assignment

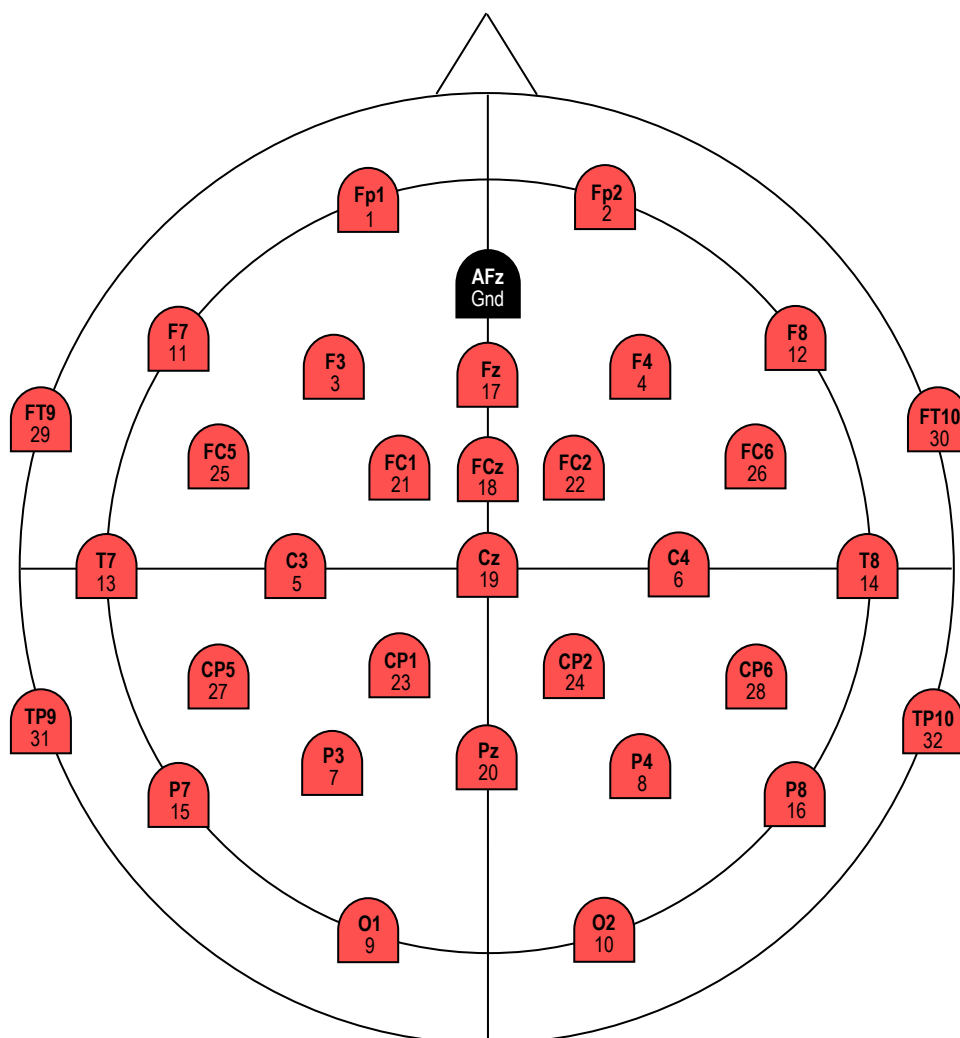

## Details

### Ordering Information

For ordering please give **Article Number, Cap Cut, and Size**  
(e.g. *BCA-32, C-Cut, 56*):

- Article Number: **BCA-32**
- Cap Cut: **C-Cut** or **A-Cut**
- Size (given in cm head circumference):
  - Adult caps: **54, 56, 58, 60, 62, 64** (average male: 58, average female: 56)
  - Children caps: **50** (3-4 years), **52** (5-10 years), **54** (11-14 years)
  - Infant caps: **34, 36** (newborn), **38, 40** (3 months), **42, 44** (7 month), **46, 48** (2 years)

The catalogue-number comprises the cap as described. For further information about accessories or consumables, please visit our website or contact our local distributor.

### Cap

Standard: Subinon Cap with integrated chin belt, white

Sizes 52 – 64 made from High Precision Fabric, Sizes 50 and smaller made from High Comfort Fabric

*Options: C-Cut or A-Cut, Size. For further variations, please contact us.*

### Electrodes

All electrodes are Multitrodes with sintered Ag/AgCl sensors. They are buttoned directly into the cap (total height 3,5 mm) or can be attached to the skin with washers (= double-sided adhesive rings).

All electrodes are number-labelled at the sensor end.

The cable colours correspond to the above figure.

The cables are attached to the cap with nylon threads. Two cable trees leave the cap plait-like behind the ears, pointing downwards, and meet before the connector-box. The length of cable trees is approx. 120 cm.

### Termination

The cable tree leads into a Connector box. From here the caps are connected to actiCHamp Plus with flat ribbon cables. Connecting this cap will require one flat ribbon cable (1x EasyCap Cat. No. *E81-22*) plus one GND connection (1x EasyCap Cat. No. *KB-B4P-B4P-20*). These cables are not included in the scope of delivery and are recommended to be ordered once per amplifier.

### Theta/Phi-Coordinates

Please find a table with Theta/Phi-Coordinates of all electrode sites at the end of this file.

Table of Coordinates for BCA-32

Channel-number	Name	Theta	Phi
1	Fp1	-90	-72
2	Fp2	90	72
3	F3	-60	-51
4	F4	60	51
5	C3	-45	0
6	C4	45	0
7	P3	-60	51
8	P4	60	-51
9	O1	-90	72
10	O2	90	-72
11	F7	-90	-36
12	F8	90	36
13	T7	-90	0
14	T8	90	0
15	P7	-90	36
16	P8	90	-36
17	Fz	45	90
18	FCz	23	90
19	Cz	0	0
20	Pz	45	-90
21	FC1	-31	-46
22	FC2	31	46
23	CP1	-31	46
24	CP2	31	-46
25	FC5	-69	-21
26	FC6	69	21
27	CP5	-69	21
28	CP6	69	-21
29	FT9	-113	-18
30	FT10	113	18
31	TP9	-113	18
32	TP10	113	-18
Gnd	Afz	67	90

These values are standardized to a Theta of 90° for the plane through Fpz, T7, T8, Oz.

The signs follow this convention:

