GACS-32-X2 - Printed 32Ch actiCAP snap for Kids (Dino)


Cap with printed prehistoric design and holders for actiCAP slim electrodes

Green Holders: Label 1-32, Hardware Channel 1 – 32

Black Holder: Label & hard-wired GND

Blue Holder: Label & hard-wired REF

Electrode Names and Number Labels

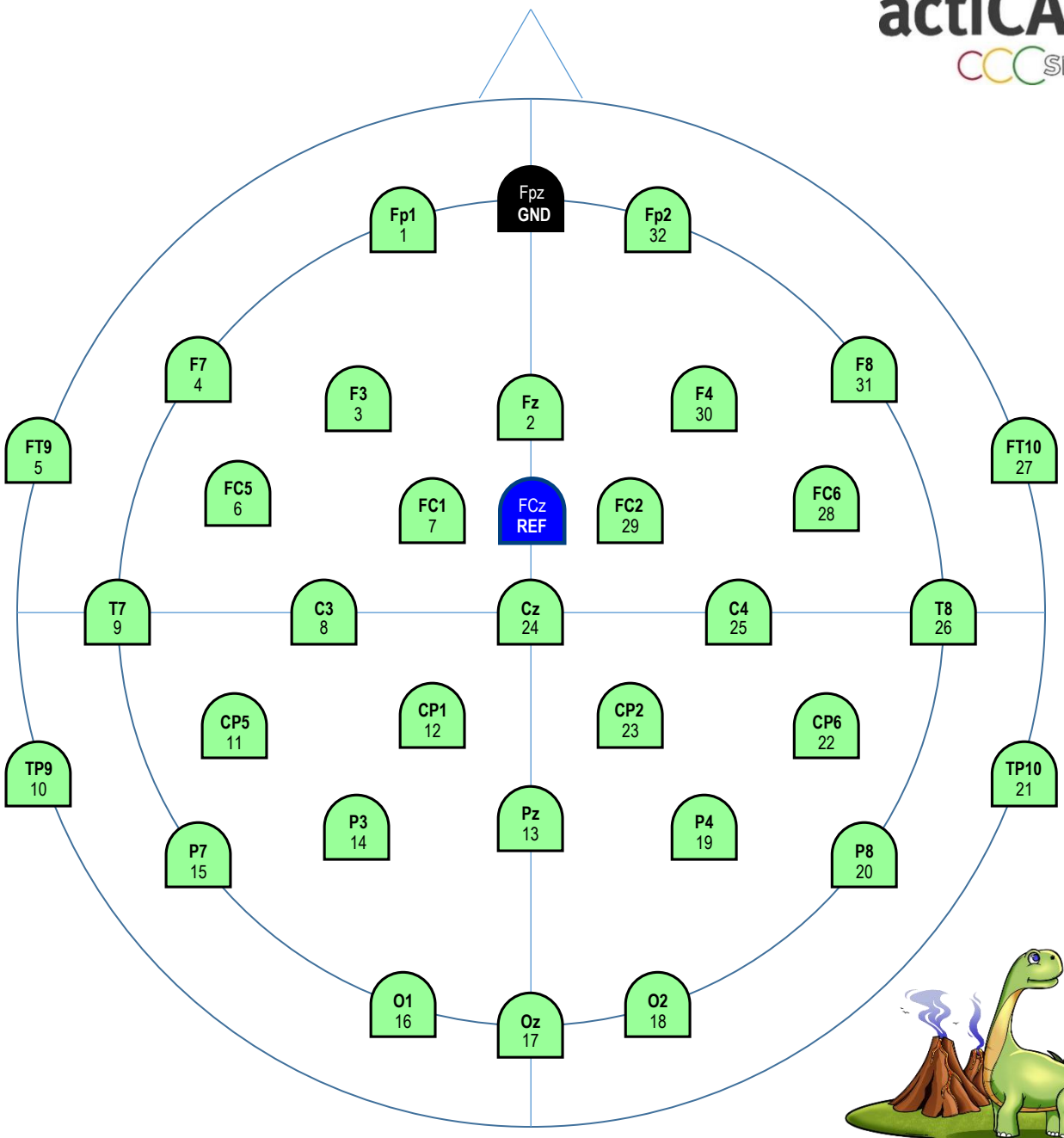



Designed by Ohmzotor (fiverr.com/Ohmzotor)

Easycap GmbH Am Anger 5, 82237 Wörthsee, Germany www.easycap.de support@easycap.de

Details



Ordering Information

For ordering please give **Layout Name, Cap Cut, and Size** *(Example: GACS-32-X2, C-Cut, 58):*

- Layout Name: ***GACS-32-X2***
- Cap Cut: ***C-Cut*** or ***A-Cut***
- Size (given in cm head circumference):
  Children caps: ***50*** (3-4 years), ***52*** (5-10 years), ***54*** (11-14 years)
  Infant caps: ***42***, ***44*** (7 month), ***46***, ***48*** (2 years)

The catalogue-number comprises the printed cap as described. For further information about accessories or consumables please visit our website or contact our local distributor.

Cap

Standard: Subinion cap with integrated chin belt, white, High Comfort Fabric, Printed Design for Kids (Dino).
Options: C-Cut or A-Cut, Size

The holders are mounted according to the figure. The positions are chosen from the 10%-System and their name labels are printed onto the cap fabric. Positions TP9/TP10 approximately meet left and right mastoids.

Apart from using the cap holders, it is also possible to use some electrodes on bare skin, e.g. for recording EOG. To do so, the electrodes can be attached to the skin with double-sided adhesive rings (washers). The grip of the washer is better, when both, skin and electrode, are degreased with alcohol prior to application.

Electrode Assignment

Plug the GND-Electrode into the black holder at Fpz.
Plug the electrodes 1-32 of the first bundle into the green holders (labelled 1-32), for amplifier channels 1-32.

If your hardware requires a REF electrode (e.g. actiCAP control box), plug REF into the blue holder at FCz.

If your hardware uses NO REF electrode (e.g. actiCHamp), leave the blue holder at FCz unused.

The velcro straps attached to the cap help you to organize and fixate the cable wires.

Cap Layout download

The QR code and link below take you to the download section of our website, where you can search for and download this cap layout document including Theta/Phi coordinates:
https://easycap.de/cap-layout-download/